

FULLY ELECTRIC
UNIVERSAL OPERATING TABLE

ST. FRANCIS has many years of extensive experience in manufacturing operating tables and operating lamps, the complete line refer to our website www.st-francis.com.tw



ST.FRANCIS Medical Equipment Co., Ltd.

P.O. Box: 129 Shien Chung, 24299 Taipei Shien, Taiwan. Office: 4F-3, No.20, Lane 609,
Sec.5, Chung shien Rd., Sanchung City, 24159 Taipei Shien, Taiwan.

Tel: 886-2-22782889.22783509.22781618 Fax: 886-2-22783497.22782799

E-mail: sthong@ms18.hinet.net


www.st-francis.com.tw



CE ISO9001. ISO13485
FDA, 510(k), GMP
CERTIFIED FACTORY

OT-2300

Fully electric operating table use for general surgery as Gynecology, Urology operation, Nephrectomy surgery and Neurosurgery procedures... is designed for your best surgeon team's requirement.

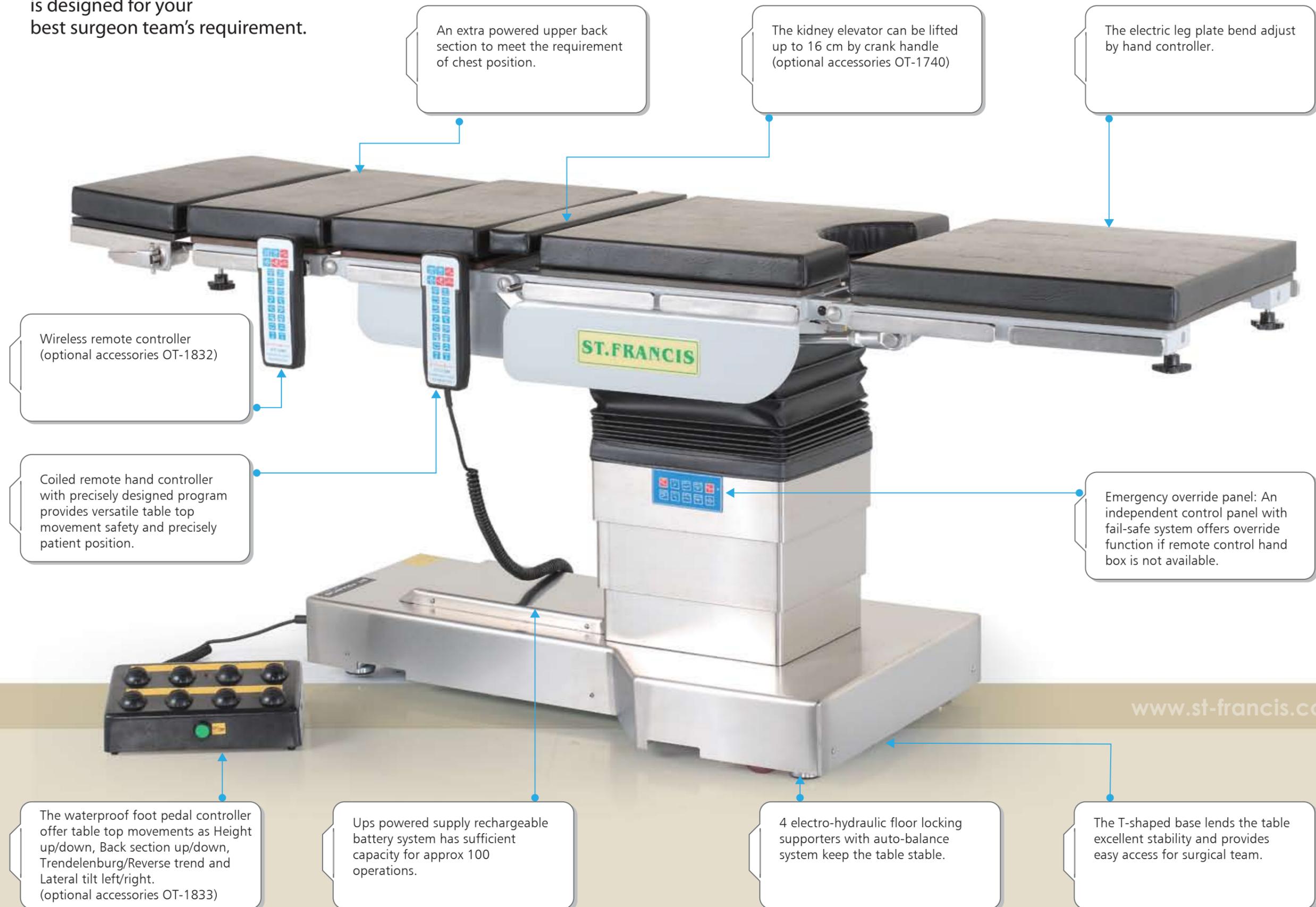

An extra powered upper back section to meet the requirement of chest position.

The kidney elevator can be lifted up to 16 cm by crank handle (optional accessories OT-1740)

The electric leg plate bend adjust by hand controller.

Wireless remote controller (optional accessories OT-1832)

Coiled remote hand controller with precisely designed program provides versatile table top movement safety and precisely patient position.

Emergency override panel: An independent control panel with fail-safe system offers override function if remote control hand box is not available.

The waterproof foot pedal controller offer table top movements as Height up/down, Back section up/down, Trendelenburg/Reverse trend and Lateral tilt left/right. (optional accessories OT-1833)

Ups powered supply rechargeable battery system has sufficient capacity for approx 100 operations.

4 electro-hydraulic floor locking supporters with auto-balance system keep the table stable.

The T-shaped base lends the table excellent stability and provides easy access for surgical team.



Vertical position for C arm x-ray



Nephrectomy position
(Optional accessories recommended: Head cushion OT-1723, Raised armrest OT-1728, Kidney Elevator OT-1740)



Thoracic position
(Optional accessories recommended: Kidney Elevator OT-1740)

Rectal position
(Optional accessories recommended: Cylindrical support OT-1726, Foot support OT-1724)



Ophthalmic position
(Optional accessories recommended: Wrist support OT-1733)



Neurosurgery position
(Optional accessories recommended: U-shaped headrest OT-1720)





Sitting position
(Optional accessories recommended: Sitting head support OT-1732, Raised armrest OT-1728)



Gynecology, Urology position
(Optional accessories recommended: Knee crutches OT-1824, Funnel and waste receptacle and corrugated tubing OT-1718)



FULLY ELECTRIC AND CONVENTIONAL MANUAL OPERATING TABLES

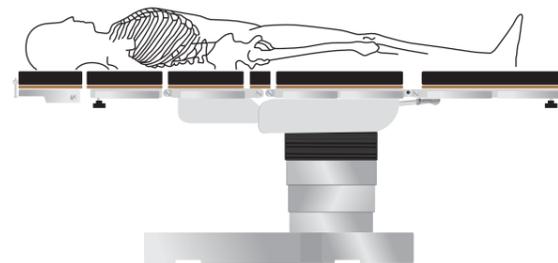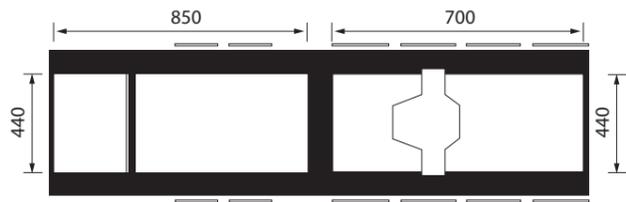


(Normal Position)

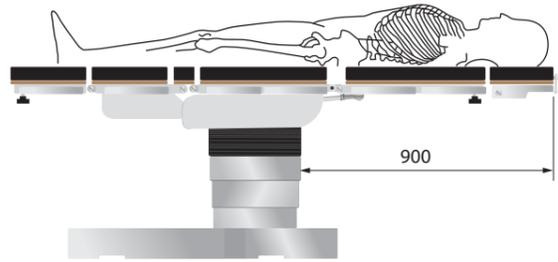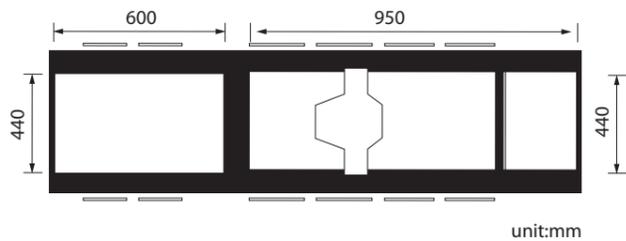

(Reverse Position / Head plate join in leg plate)

▶▶▶ OT-2300
STANDARD ACCESSORIES

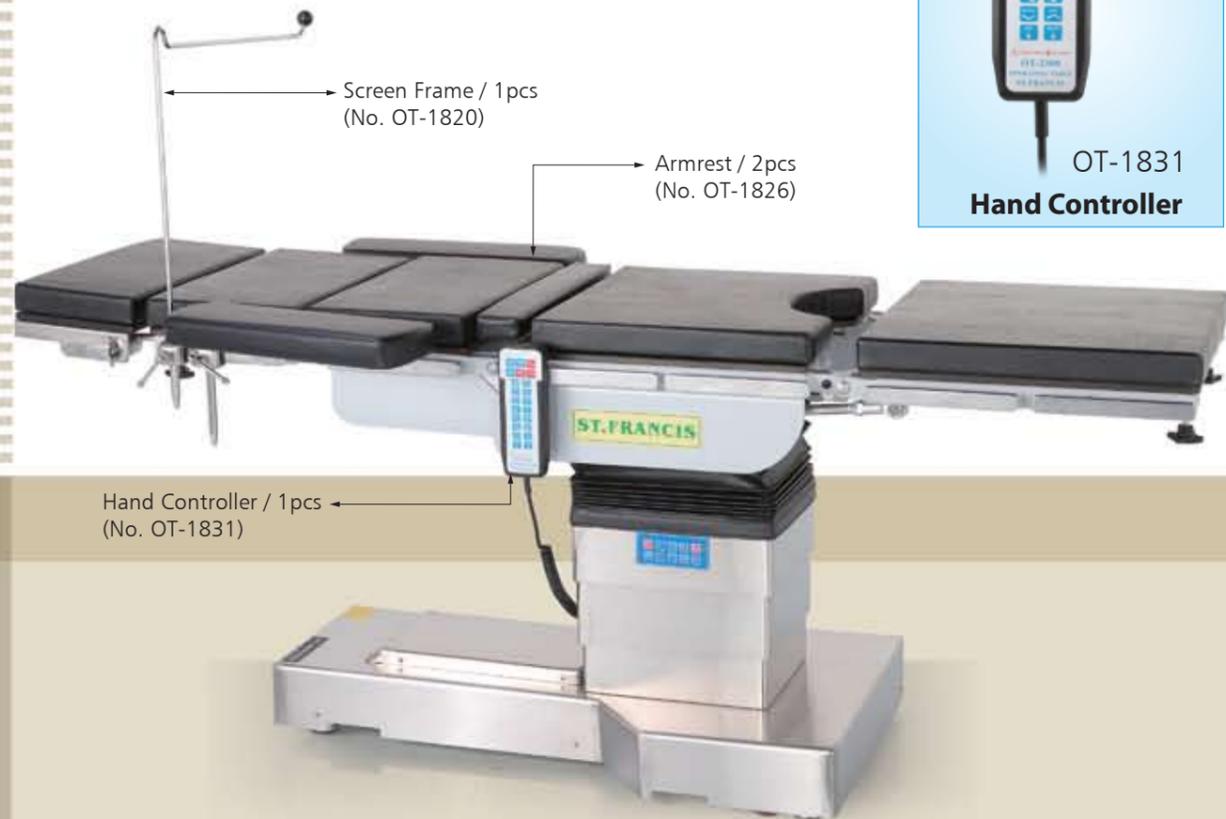

Screen Frame / 1pcs
(No. OT-1820)

Armrest / 2pcs
(No. OT-1826)

Hand Controller / 1pcs
(No. OT-1831)



OT-1831
Hand Controller

Model Number	OT-2300	
Descriptions		
Table Dimension/mm		
Table Top Elevation		
Trendelenburg/ Reverse Trend		
Lateral Tilt (Left/ Right)		
Lower Back Section (Up/Down)		
Upper Back Section (Up/Down)		
Head Section (Up/Down)		
Footplate (Up/Down)		
Flex/Reverse Flex		
Emergency Override Panel	Standard Accessories	
Ups Battery Device	Standard Accessories	
Elec. AC Input	110V ~ 230V 50/60Hz	
Remote Controller Power	DC 12V	
Table Weight (Approx)	285	
Mattress Thick	5.5 CM	
Patient Weight Capacity	250 KGS	