



OT-2200 . OT-2100

FULLY ELECTRIC UNIVERSAL OPERATING TABLE

The latest multi-function operating table is providing general surgery as:
Urology, Gynaecological, E.N.T., Ophthalmic, Rectal, Neurosurgical and Orthopaedic procedures.
CE, ISO13485, FDA 510(k), GMP certified factory.

Advance Technology in Surgical Table

- Table top consists of five sections covered with detachable electric conductive fiber mattress to provide high performance flexibility. The head and leg boards can be removed manually.
- Touch-Screen coiled remote hand controller with precisely designed program provides versatile table top movement safety and precisely patient position.
- 4 castors can roll both longitudinally and laterally.
- The T-shaped base lends the table excellent stability and provides easy access for surgical team.
- 4 electro-hydraulic floor locking supporters with auto-balance system keep the table stable.
- Table top longitudinal slides up to 35cm to get the maximum radiolucent field for free access C-ARM without moving the patient or extra extensions on table top.



Standard Accessories


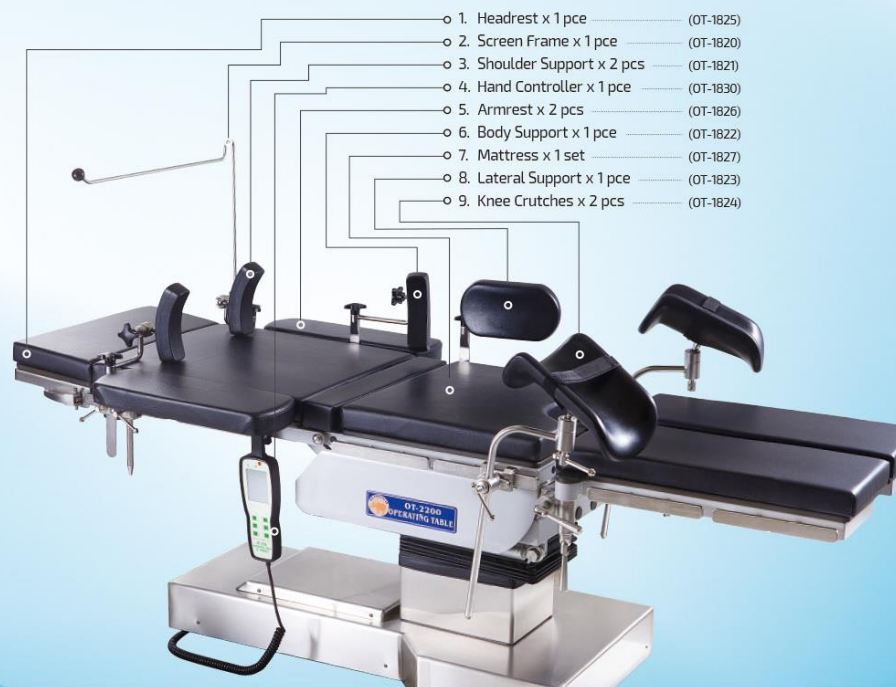


Table top sideways sliding 15cm provides a perfect position for X-ray radiography and surgery procedures.



Two separate boards can be adjusted up, down or split tilt angle position.

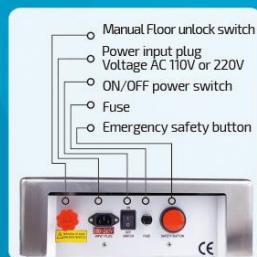

UPS powered supply rechargeable battery system has sufficient capacity for approx 100 operation. (Standard device for OT-2200)



X-ray photography table top offers whole body image amplification coverage.



Emergency override panel: An independent control panel with fail-safe system offers override function if remote control hand box is not available. (Standard device for OT-2200)

Advance Technology in Surgical Table

OT-2200/2100

FULLY ELECTRIC UNIVERSAL OPERATING TABLE



Table Dimension Descriptions

Descriptions	Model	OT-2200	OT-2100
Table Dimension (Length/Width)		2125 mm / 600 mm	
Table Top Elevation		740~1040 mm	
Trendelenburg/Reverse Trend		20°/20°	
Lateral Tilt (Left/Right)		20°/20°	
Back Section (Up/Down)		80°/30°	
Head Section (Up/Down)		25°/90°	
Footplate (Up/Down/Split Tilt)		20°/90° (Manual)	
Table Top Longitudinal Slide		350 mm	
Table Top Sidewise Slide (Left/Right)		150 mm	N.A
Flex/Reverse Flex		210°/90° (Auto-Setting)	210°/90° (Manual)
Emergency Override Panel		Standard	
Ups Battery Device		Standard	Optional
Elec. AC Input		110~230 V, 50/60 Hz	
Remote Controller Power		DC 12V	
Table Weight (Approx)		285 kgs	
Mattress Thick		55 mm	
Patient Weight Capacity		250 kgs	
Safety Standards		CE, FDA, GMP	

The specification and design to change without notice.



ST. FRANCIS MEDICAL EQUIPMENT CO., LTD.

4F-3, No. 20, Lane 609, Sec.5, Chung Shien Rd.,
Sanchung District, New Taipei City, Taiwan
Tel: 886-2-22782889, 22783509
Fax: 886-2-22783497, 22782799
sthong@ms18.hinet.net
www.st-francis.com.tw



CE, ISO13485, FDA 510(k), GMP certified factory

The specification and design to change without notice.
Printed in Taiwan, May, 2018 All Right reserved.
St. Francis Catalogue Ref. OT2200-E2



OT-2200 . OT-2100

FULLY ELECTRIC UNIVERSAL OPERATING TABLE

The latest multi-function operating table is providing general surgery as:
Urology, Gynaecological, ENT, Ophthalmic, Rectal, Neurosurgical and Orthopaedic procedures.
CE, ISO13485, FDA 510(k), GMP certified factory.

